

A black and white photograph of a woman with short hair, wearing a long-sleeved top and dark pants, standing next to a curved stairlift. She is smiling and has her hand on the backrest of the chair. The stairlift is white and is positioned on a carpeted staircase with a wooden handrail. The background shows the staircase curving upwards.

Stairfriend 23

Curved Stairlift

Owner's
Manual

 **savaria**

IMPORTANT – READ FIRST

Ensure that only an authorized Savaria mechanic installs and services your Lift. Under no circumstances is anyone other than a mechanic with Savaria training and authorization to install, adjust, service or modify any mechanical or electrical device on this equipment.

Failure to follow this warning can result in safety system compromises or defeat; this can result in serious injury or death. Savaria accepts no liability for property damage, warranty claims or personal injury, including death, in this circumstance.

Passenger safety is the result of countless details in the equipment design, manufacture, and installation. After installation, reliable operation and continual safe operation requires regular inspection and maintenance as outlined in this manual, or more frequently where usage, environment, or local jurisdiction requires. As the Owner, you are responsible for ensuring that regular inspections and maintenance occur in a timely manner as per the manufacturer's recommendations and maintenance schedule.

Refer to this owner's manual for unit description, setup, operating instructions, and a maintenance schedule for your Lift. Upon completion of installation, the Savaria mechanic must provide you with the information for the box below and ensure it is recorded in this owner's manual. In addition, either you or the Savaria mechanic must keep any service and/or maintenance records in the Maintenance Record section of this owner's manual.

At all times, your Lift must be maintained in compliance with the applicable standard to which it has been built and installed.

United States: A18.1–Safety Standard for Platform Lifts and Stairway Chairlift

Canada: CSA B355/613–Platform lifts and stair lifts for barrier-free access

Europe: EN-81-40–Stairlifts and inclined lifting platform intended for persons with impaired mobility

Europe: 2006/42/EC–Machinery Directive

Depending on your country or region, the local authority code variant may apply. To ensure safe operation of your Lift, pay careful attention to the important notes, and follow the instructions and any additional safety regulations in this manual.

FOR OWNER'S RECORDS

Customer Name: _____

Installing Dealer: _____

Dealer's Telephone Number: _____

Date Installed: _____

Serial/Job Number: _____

CONTENTS

1. SAFETY INSTRUCTIONS	4
To Ensure Safety	4
2. DESCRIPTION	5
Seating Features	7
Safety Features	7
Footrest and Carriage Sensors	7
Safety Switches	7
Emergency Stop Button (Optional)	7
Overspeed Governor	7
Key Switch and Power Switch	8
3. SPECIFICATIONS	9
4. USING THE DEVICE	10
Before Operating the Stairlift	10
Operating the Stairlift	11
Emergency Stop	12
Emergency Manual Operation.....	12
5. ADJUSTMENTS	13
Seat Adjustment (V-Seat)	13
6. TROUBLESHOOTING	14
Diagnostic Codes	15
6. MAINTENANCE	19
Maintenance Schedule	19
Maintenance Record	23

1. SAFETY INSTRUCTIONS

To ensure safe operation of this unit, pay careful attention to the important notes below.

To Ensure Safety

- Read this manual carefully before using the equipment.
- To prevent accidents, adhere strictly to the instructions and keep clear of moving parts at all times.
- Follow instructions on all equipment labels at all times. Replace any damaged labels immediately.
- Ensure that only qualified personnel perform maintenance and service on the unit.
- When replacing parts, be sure that only genuine Savaria parts are used.
- This unit is intended for use by a mature person who understands its proper operation as set out in this manual.



WARNING

At no time should children under the age of 16 use the lift while unsupervised.


- This lift is to be used **ONLY** for the transport of a person up/down the stairway. Do **NOT** use the lift for any other purpose (such as the transport of items).
- Do **NOT** use the lift for firefighting purposes or for evacuation during a fire.
- Prior to operation, make sure that all areas in and around the lift are clear of any obstructions.
- Be very careful of others using the stairs while the lift is in operation as there is a potential risk of collision.
- Test your keys and emergency stop button (if equipped) every month.

2. DESCRIPTION

Exterior View

Figure 1 shows the exterior components of the Stairfriend 23 with V-seat.

Figure 1

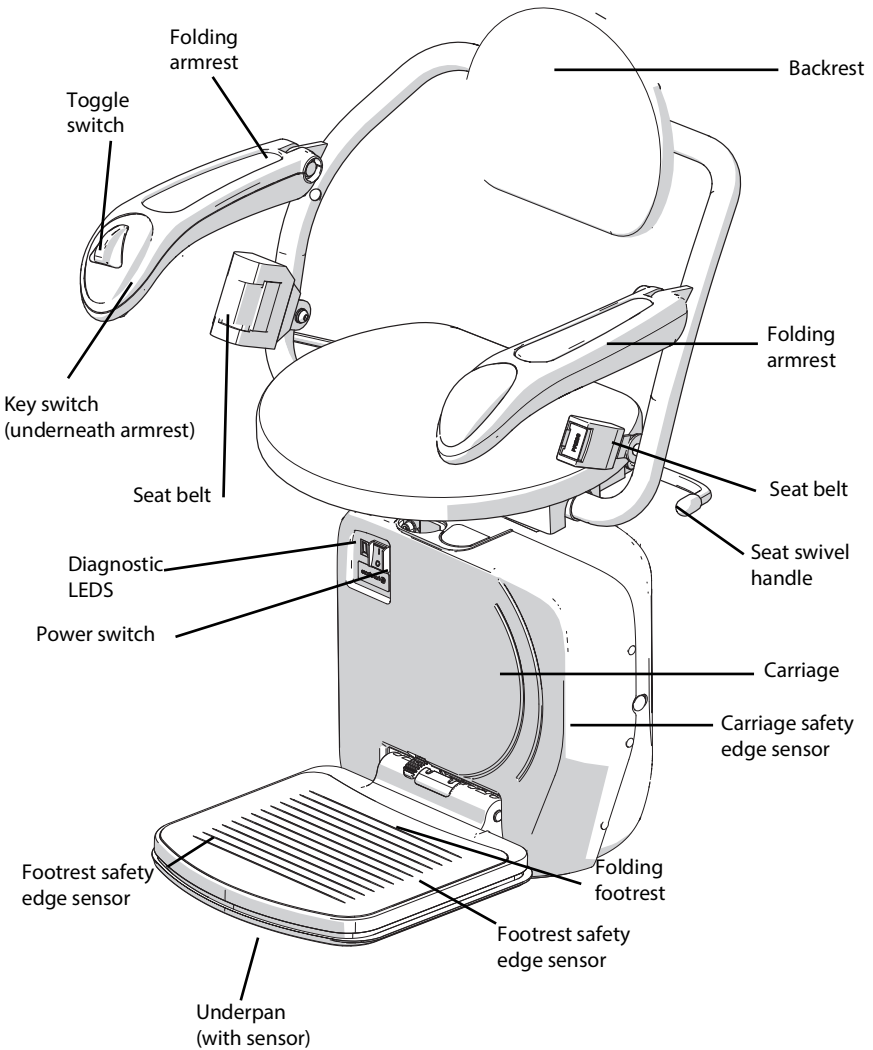

Seating Features

The seat assembly can swivel in either direction for transfer in/out of the seat. Lift up and hold the swivel handle to swivel the seat. Let go of the handle and the seat will lock in place. **For your safety, do not swivel the seat until the unit has stopped.**

When the stairlift is not in use, the seat, armrests, and footrest can be folded up to allow for more space on the stairway.

Always use the retractable seat belt to ensure that you are safely secured in the chair during travel.

Safety Features

Footrest and Carriage Sensors

There are safety edge sensors on the sides of the footrest and the carriage, and an underpan sensor on the bottom of the footrest. These sensors will stop the stairlift if an obstruction is encountered on the stairway.

If you encounter an obstruction, press the opposite side of the toggle switch to run the stairlift in the opposite direction. Remove the obstacle from the pathway and proceed in the original direction of travel.

Safety Switches

There are upper and lower safety limit switches to stop the stairlift when it reaches the top or bottom of the stairway. There are also mechanical stops at the top and bottom of the rail to stop the lift from running off the rail.

There is also a safety switch for the swivel seat that will prevent the stairlift from running if the seat is not locked in the riding position.

Emergency Stop Button (Optional)

There is an optional emergency stop button that may be located on the arm. Pressing this button will stop the unit immediately.


Overspeed Governor

There is an overspeed governor that will engage and stop the lift if the normal downward travel speed is exceeded.

Key Switch and Power Switch

The stairlift is equipped with a key switch underneath the armrest (see Figure 3). The stairlift cannot be operated unless the key is inserted and turned to the ON position. Once the key switch is turned on, flip the red ON/OFF power switch located on the front top left of the carriage to the ON position.

Figure 3

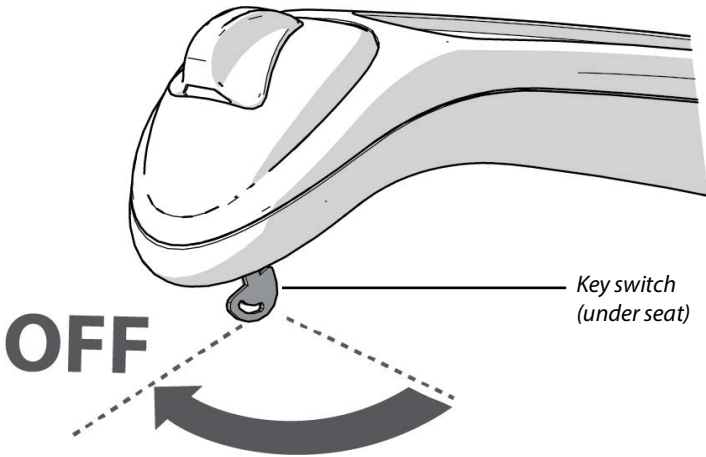

3. SPECIFICATIONS

Stairfriend 23 Specifications

Applications	Residential, indoor, curved stairs Commercial - varies by region; contact your authorized Savaria dealer for details
Load capacity	300 lb (136 kg), Heavy duty, 350 lb (160 kg)
Maximum travel	49 ft (14.9 m)
Travel speed	Approximately 20 ft/min (0.1 m/s)
Daily cycle	Normal: 16 Maximum starts in 1 hour on standard installation: 10
Range of incline	0° to 56°
Drive system	Rack and pinion
Operator controls	Constant-pressure type (toggle switch and remote control call/send units)
Seat belt	Retractable
Motor	1/2 HP, 24 VDC
Charging	Charges on standard household power (110 VAC) at landings
Fold-up width	16.25" (41.2 cm)
Fold-out width	27" (68.6 cm)
Rail	Steel
Emergency operation	Full up and down travel on power interruption (full battery backup)
Options	Savaria Link remote monitoring
Warranty	36 months parts, ask dealer for details

4. USING THE DEVICE

Before Operating the Stairlift

- 1 Ensure the stairway is clear of any obstructions and the stairlift charger is plugged in at the wall.
- 2 Insert the key and turn it a quarter turn clockwise. Flip the ON/OFF power switch to the ON position.
- 3 Check that the seat is locked in the riding position for up/down travel.

NOTE

If the chair has not been used for over an hour, press the remote button twice to wake it from sleep mode. If the chair has not been used for over 24 hours, run the chair from the arm toggle switch for a second to wake it from sleep mode. You can then use the remotes for normal operation.

Operating the Stairlift

- 1** Use one of the remote control CALL/SEND units to move the stairlift to the required end of the rail. Be sure to point the remote at the lift when pressing the button.
- 2** If necessary, fold down the seat, armrests, and footrest.
- 3** Lift up the seat swivel handle, rotate the seat to a comfortable boarding position, and then release the swivel handle to lock the seat in place.
- 4** Position yourself in the seat and fasten the seat belt.
- 5** Lift up and hold the seat swivel handle and rotate the seat so you are parallel to the rail (with your back against the stairway wall). Release the swivel handle to lock the seat in the riding position.
- 6** Place your feet securely on the footrest and your arms on the armrests.
- 7** Press and hold one end of the toggle switch to move the stairlift in the desired direction.
- 8** The stairlift will stop automatically at the upper or lower end of the rail, or if the toggle switch is released.
- 9** Once the stairlift has come to a complete stop, swivel the seat as required. When the seat is locked in place, unfasten the seat belt and get out of the chair.
- 10** To allow for more room on the stairway when the stairlift is not in use, swivel the seat back to the riding position, and then fold up the seat, armrests, and footrest.

Emergency Stop (Optional)

In the event of an emergency, you can press the optional red Emergency Stop button located on the arm of the chair. The stairlift will stop immediately. Once the emergency is dealt with, release the button to restore normal operation.

Emergency Manual Operation


The emergency manual lowering tool (hand winder) allows you to operate the lift manually.

NOTE

The batteries must be connected to release the motor brake and then you can use the manual hand crank.

- 1 Turn off power at the stairlift and unplug the power cord from the wall.
- 2 Remove the rubber cover from the hole on the top front of the carriage.
- 3 Fit the shaft of the manual lowering hand winder onto the motor input shaft.
- 4 Rotate the tool **clockwise** for the unit to travel **down** the stairs and **counter-clockwise** for the unit to travel **up** the stairs.
- 5 When complete, remove the tool and insert the rubber plug back into the hole.

Figure 4

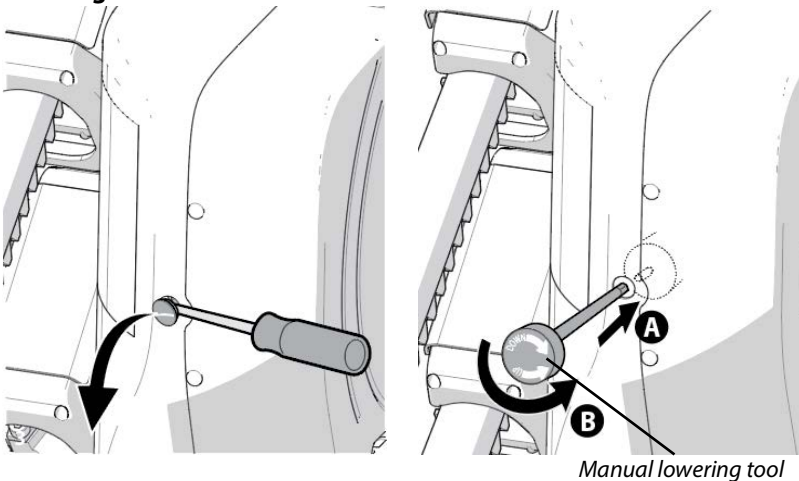



5. ADJUSTMENTS

Seat Adjustment (V-Seat)

If the seat will not stay folded up when the stairlift is not in use, perform the adjustment procedure below.

- 1 Lift up the seat and loosen the two lock-nuts as shown below.
- 2 Using an Allen key, tighten the set-screw in each lock-nut as shown below.
- 3 After tightening the set-screws, move the seat to ensure it will remain in the folded position.
- 4 Tighten the two lock-nuts.

Figure 5

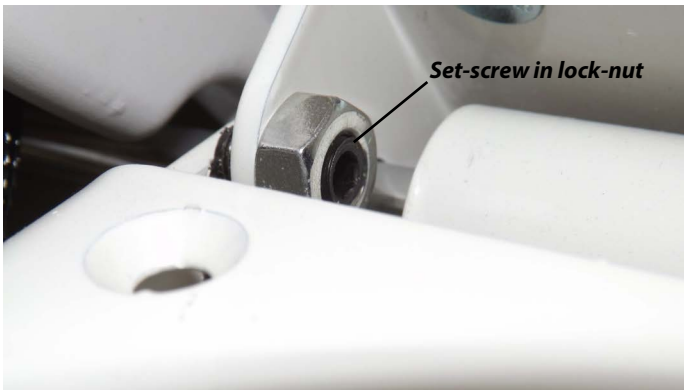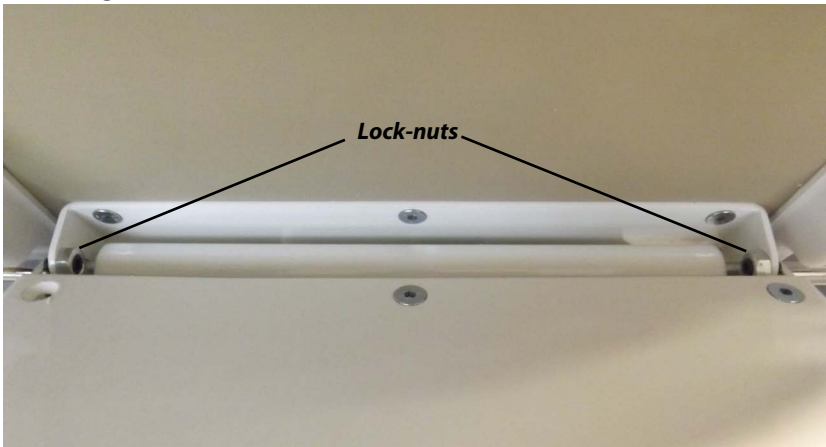

6. TROUBLESHOOTING










WARNING









Only qualified technicians are authorized to perform repairs and maintenance. DO NOT manipulate, modify or remove any safety feature of the lift. If any issues persist, please contact your authorized Savaria Dealer









Problem	Action
Stairlift will not run and no lights are on	<p>Make sure the stairlift charger is plugged in at the wall.</p> <p>Make sure the key switch is turned ON and the power switch is ON.</p> <p>Check the residence main circuit breaker (or fuse) panel for a tripped breaker tripped or blown fuse?</p>
Stairlift does not operate and all of the lights are on	<p>-Make sure the seat is swiveled and locked in the riding position.</p> <p>Check around the footrest and carriage for any obstructions.</p>
Stairlift is beeping	<p>If the stairlift beeps for approximately 20 seconds while not moving, this means the stairlift is not charging.</p> <p>Check that the stairlift cord is plugged in to the wall outlet.</p> <p>Make sure the stairlift is either all the way up or all the way down at the charging contact. It will beep continuously if it is not at a charging contact.</p>
Stairlift stops and won't restart	<p>f the stairlift runs to either end of the stairway and then stops and won't restart, this means that the battery charge is low.</p> <p>Charge the stairlift.</p>



Diagnostic Codes

Figure 6

	No Display
	OSG Input Low
	Off charge.
	Top (LH) or Bottom (RH) track limit activated.
	Safety Up
	Under Pan
	Safety Down

	Low battery voltage. Requires charge.
	Switch active at Power On
	Toggle Up or Down
	IR address fail.
	Power relay closed and should be open (welded)
	Power relay Open and should be closed
	Final limit low
	Brake open and should be closed

	Power relay closed early
	Ultimate limit
	Over current
	Slow speed
	Memory problem.
	No charge current
	FET Fail
	Seat swivel

	Key Switch
	Power supply fault when battery on float

6. MAINTENANCE

Stairfriend 23 is subject to wear and tear from use. You must perform the checks and actions in the following tables to ensure safety and proper operation.

IMPORTANT

Savaria products are only to be installed, adjusted, serviced, or maintained by Savaria licensed dealers and technicians. Your Savaria product will have the warranty voided if a non-Savaria approved technician performs work on the Savaria product.

For units with high or excessive daily cycles, outdoor use, or use in harsh environments, inspection and maintenance should be conducted more frequently to ensure optimal performance.

Maintenance Schedule

Verification	Frequency (Commercial)	Frequency (Residential)
Clean off the gear rack with a dry cloth and then apply the clear, multi-purpose grease	Every 6 months	Normal: Every Year Heavy: Every Year Excessive: Every 6 months
Check that the swivel seat switch is working. Ensure the unit will not travel until the swivel seat is locked in the riding position.	Every 6 months	Normal: Every Year Heavy: Every Year Excessive: Every 6 months
Check that the footrest and carriage sensors are working. Ensure the unit stops when it contacts an obstacle.	Every 6 months	Normal: Every Year Heavy: Every Year Excessive: Every 6 months

<p>Check that the emergency stop button is working (if equipped). Ensure the unit stops immediately when the button is pressed.</p>	<p>Every 6 months</p>	<p>Normal: Every Year Heavy: Every Year Excessive: Every 6 months</p>
<p>Check the operation of the limit switches. Ensure the unit stops at the required position at the top and bottom of the stairway.</p>	<p>Every 6 months</p>	<p>Normal: Every Year Heavy: Every Year Excessive: Every 6 months</p>
<p>Check the operation of the overspeed governor</p>	<p>Every 6 months</p>	<p>Normal: Every Year Heavy: Every Year Excessive: Every 6 months</p>
<p>Lubricate the trolley rollers with silicone grease. Refer to the Installation Guide Maintenance section for the procedure.</p>	<p>As required (if squeaking occurs)</p>	<p>As required (if squeaking occurs)</p>
<p>Lubricate the seat swivel system using wheel bearing grease.</p>	<p>As required</p>	<p>As required</p>

THIS PAGE WAS INTENTIONALLY LEFT BLANK

THIS PAGE WAS INTENTIONALLY LEFT BLANK

Stairfriend 23

Curved Stairlift

Owner's Manual

For service or questions about this product,
please contact your installing dealer.

Dealer Name: _____

Dealer Phone: _____

Authorized Savaria Dealer

Savaria Concord Lifts, Inc.
2 Walker Drive Brampton ON L6T 5E1 Canada



savaria.com